

With over thirty years of experience, Carilex has established itself as a global leader in the wound care industry by developing state-of-the-art technology in the prevention and treatment of pressure ulcers. Established relationships with medical institutions worldwide allow us to gain critical insight into the needs of our patients and their clinicians. Our expertise in research and development; manufacturing, and quality control and assurance allow us to provide unparalleled products in wound care management in the form of comprehensive product lines proven to prevent and treat of pressure sores and wound. Our strength is in our people and their dedication to and belief in the excellence of the Carilex product. It is evident in our continuous and rigorous devotion to providing the finest medical devices.

Together with commitment and constant innovation, Carilex is driven forward by our motto—CARING MAKES THE DIFFERENCE. All of us at Carilex believe paying extra care in everything we do is the key to make real differences.

We take pride and find joy knowing that our patients are in good hands when they are using Carilex products.

EUROPE, MIDDLE EAST & AFRICA
Carilex Medical GmbH
Turmbergstr, 9A, D-76227 Karlsruhe, Germany

TEL : +49 (0) 721 16 19 835
emea@carilexmedical.com

UNITED KINGDOM
Carilex Medical UK
22-23 Mountney Bridge Business Park, Westham, East Sussex,
England, BN24 5NH

TEL : +44 (0) 1323 761900
uk@carilexmedical.com

ASIA PACIFIC
Carilex Medical, Inc.
No.77, Keji 1st Rd., Guishan Township, Taoyuan County (333),
Taiwan (R.O.C)

TEL : +886 (0) 3 328 7882
asia@carilexmedical.com

ANZ
Carilex Medical ANZ Pty Ltd
Unit 23/10 Victoria Ave Castle Hill NSW 2154

TEL : +61 1300 732 004
anz@carilexmedical.com

AMERICAS
Carilex Medical Americas
10163 Patriot Drive, Suite 2 Baton Rouge, LA 70816

TEL : +1 855 721 7233
americas@carilexmedical.com

carilex[®]

CozinyPlus

Gentle healing for the most precious person in the world.



Gentle healing for the most precious person in the world.

Here at Carilex, our motto is "caring makes the difference". The care for our patients has always been the primary drive for us to strive for excellence. We pride ourselves in creating products that take care of our patients whenever they needed. Carilex CozinyPlus offers wound care therapy designed specifically to take care of children who are at risk of developing pressure ulcers, or already enduring distress and discomfort from pressure ulcers. CozinyPlus provides multiple mattress sizes to choose from, it is compatible with various body types and an array of environments. Details were considered thoroughly in the design of our mattresses ensuring the most gentle healing environment for sensitive skin. With teams of research and design experts working tirelessly throughout the years, CozinyPlus is a labour of love. Carilex cares deeply in what we build. We care about making a real difference in providing attentive products to those who are in need.

The CozinyPlus Power Unit with IPS Technology.

The IPS Technology, Intelligent Pressure Sensing Technology, responds to patient movements on mattress by automatically adjusting internal mattress pressure, which allows CozinyPlus to regulate airflow and continuously provide total envelopment for our patients. Carilex is the first brand in the industry to master IPS. This gives us more experience and vast knowledge in the technology that is years ahead of our competitors.

1. Multiple compatibility

Compatible with multiple mattress sizes for various body types.

2. Intelligent pressure sensing

Responds to patient movements on mattress by automatically adjusting internal, mattress pressure.

3. Multiple weight setting

Selectable weight settings to ensure optimal pressure relief for all body types.

4. Dynamic / static therapy modes

Provide both dynamic and static functions to meet individual patient's need.

5. Multiple cycle times

Various cycle time selection for alternating pressure therapy.

6. Battery backup upto 8 hours


Automatic power detection sensor activates battery backup in the event of main power failure and allows treatment to continue uninterrupted during patient transportation and examination.



The Coziny Mattresses

CozinyPlus offers multiple types of mattress designed with extra attention to details, providing remarkable comfort. We gave extra considerations in eliminating sharp or protruding components from mattresses; ensuring that we deliver the most protective wound care prevention and therapy for the most sensitive skin.

The Carilex Coziny 300 Mattress for Baby Cot

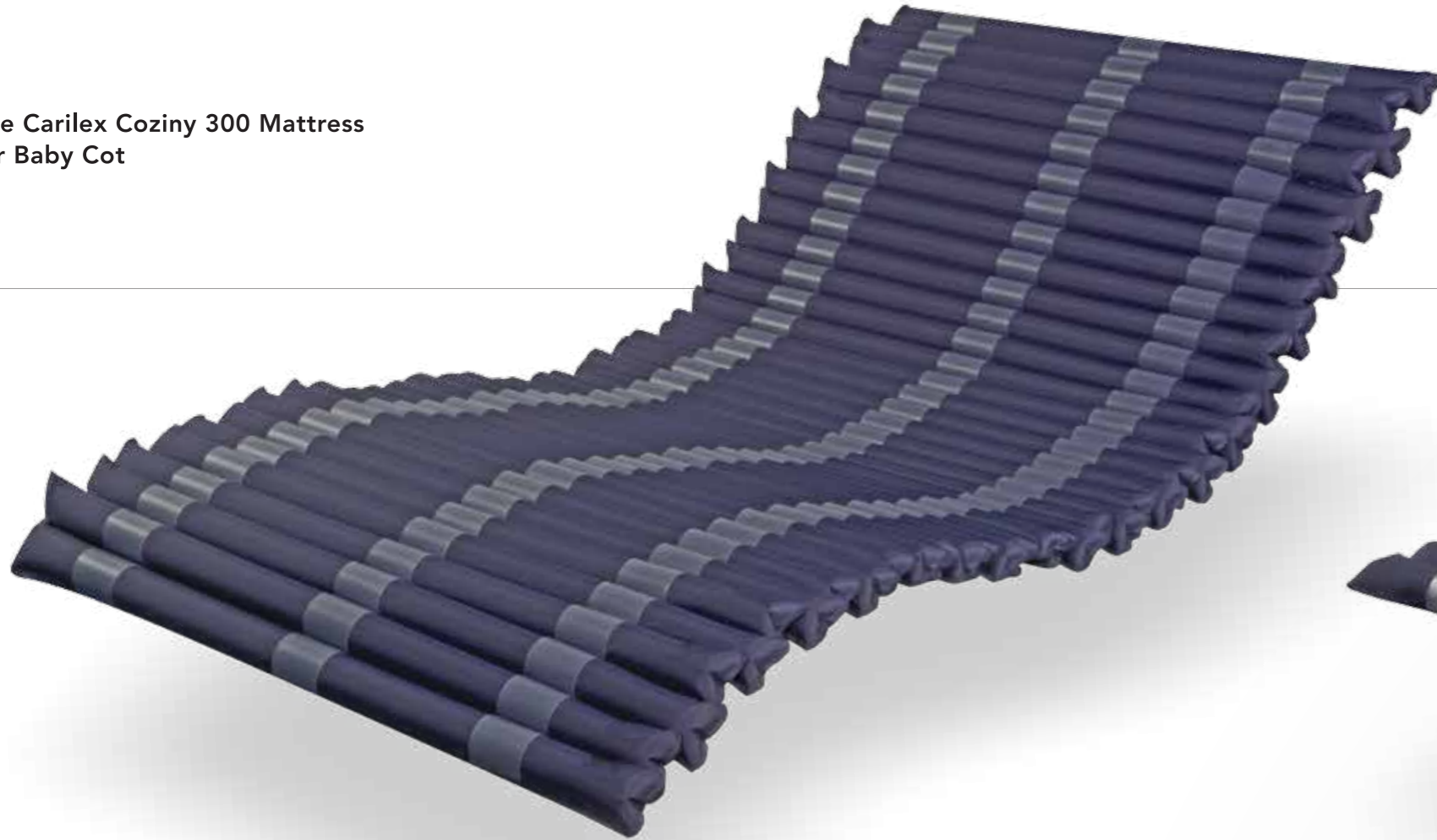

1. Three connecting hose outlets

Different positions for connecting hoses, accommodating different condition for baby cots seamlessly.


2. 2.8cm foam base

Integrated for extra protection in case of power outage

3. Weight capacity

For patient up to 25kgs

The Carilex Coziny 200 Mattress for Baby Crib

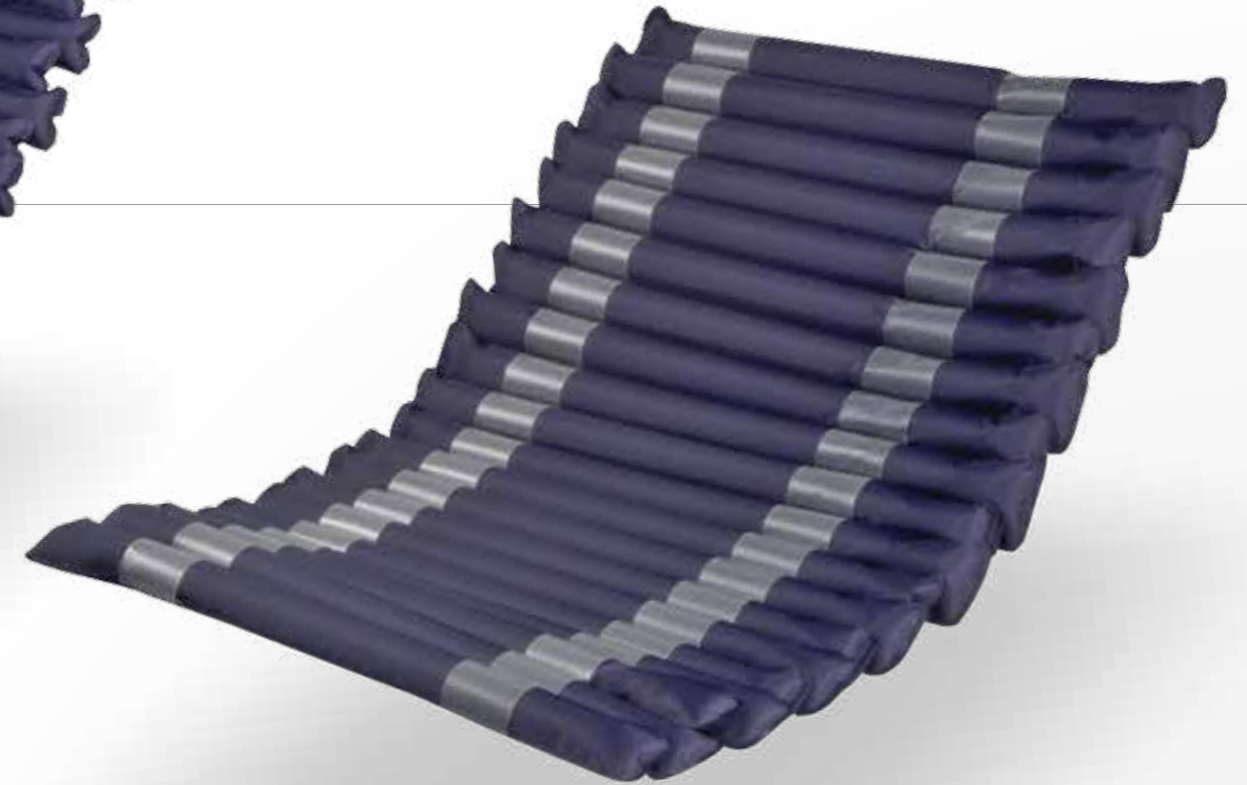

1. Three connecting hose outlets

Different positions for connecting hoses, accommodating different condition for baby cribs seamlessly.

2. 1.5cm foam base

Integrated for extra protection in case of power outage.

3. Weight capacity

For patient up to 10kgs.

The Carilex Coziny 100 Mattress for Incubator

Coziny 100 mattresses are made entirely with Extra Soft 100% TPU material to ensure baby's safety and comfort. The Extra Soft 100% TPU fabric air inlets also protects baby from sharp edges and protruding objects.

1. Extra soft 100% TPU material

Ensuring patients safety from sharp or protruding objects while providing optimal comfort.

2. Special temperature distributing air channels

Allowing the heat from incubator to evenly keeping the baby warm throughout.

3. Hidden zipper

With protection pad for baby's safety and better comfort.

4. Weight capacity

For patient up to 3.4kgs.



Other Features



LED



Panel Lock



Alarm Mute



Whisper Quiet



Visual and audible indicator



Battery level LED



Quick Connect Coupling



Modular air cell construction

Technical Data

- Power Unit Model : _____ REF : SR321
- Power Input : _____ 100-240Vac 50-60Hz 1.5A
- Power Consumption : _____ 27W_Adaptor + Pump; 9W_Pump only
- Weight : _____ Power Unit : 0.92kg / 2.03lb
 Mattress Set :
 Coziny 100, 0.5kg / 2.87lb
 Coziny 200, 1.5kg / 5.51lb
 Coziny 300, 3.5kg / 9.09lb
- Dimension : _____ Power Unit : 24.5 x 15.5 x 6.1 cm
 Mattress Set :
 Coziny 100 : 25" x 14" x 1" (65 x 35 x 3 cm)
 Coziny 200 : 30" x 20" x 3" (75 x 50 x 7 cm)
 Coziny 300 : 55" x 28" x 3" (140 x 70 x 8 cm)
 other size upon request
- Agency Approval : _____ SGS-Q Mark UL60101-1
- Operating Condition : _____ Temperature : 5°C to 40°C
 Humidity : 15% to 60%
- Transportation & Storage Condition : _____ Temperature : -25°C to +70°C
 Humidity : 0% to 93%
- Weight Capacity : _____ Coziny 100 : 3.5kg
 Coziny 200 : 10kg
 Coziny 300 : 25kg
- This Power Unit is not AP/APG protected