

PERFECT FOR PATIENTS
PERFECT FOR CAREGIVERS

Rhythm·Turn

Rhythm·Turn offers continual lateral rotation therapy to help mobilize secretions and other bodily fluids during pulmonary treatment related to immobility.



Various top covers are designed and available to satisfy your individual requirement in terms of therapy, prevention and patient comfort.

Master Control Unit

- Dimension: (L)17.7"x (W)6.8"x (H)10.8"
- Weight: 13.8 lbs
- Rated Voltage: 110-120 volts
- Max Current: 5A
- Rated Frequency: 60Hz
- Cycle Time: 3-95 mins
- 3 in 1 alternation
- Min/Max Pressure: 10-61 mmHg +/- 3-4 mmHg
- Protection Type: Class I, Type BF
- Power outage visual and audio alarm
- Static Function
- Auto Firm for fast mattress inflation
- Keypad Lock-out Function

Mattress Replacement

- Model No: Rhythm·Turn
- Dimension: (L)80"x (W)36"x (H)10"
- Weight: 37.4 lbs
- Weight Capacity: 600 lbs
- Cell Material: Blue Nylon with PU Coating
- Cell Number: 18 cells
- Cover Material: Blue Nylon with PU coating finish
- Base Material: Black Polyester with PVC Coating
- CPR Quick Release System



Rhythm·Turn



Rhythm·Turn

provides you, the carer, with peace of mind
and your patient with world class care & comfort

Perfect for
patients
and
caregivers

- High vapor permeable and quilted nylon cover

- Low Air Loss Therapy provides extra comfort for moisture and temperature control

- The Rhythm·Turn replacement mattress can accommodate patient weights up to 600 lbs

- 20°C and 40°C lateral rotation with side perimeter and static therapy




Rhythm·Turn
is your partner
in protection
and treatment
of pressure ulcer



CPR quick
release system deflates
the system
within 20 seconds.



Cable management
keeps
electrical cables
off the floor

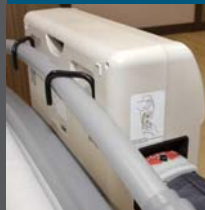

Hanging bracket
fits on
all type of beds



Low slip base cover
is a plus to keep
the mattress
in place
on the bed



Effective
micro low air loss
for moisture
management



Nylon woven
fabric with
PU coating finish